

Power wheelchairs

M1



permobil

Permobil tough

A true Permobil quality solution for those users who require a more traditional rehab style seating system with limited to no power seat functions.

Designed with adaptable seating and positioning needs in mind, the M1 offers all the flexibility without the compromise.

1" inch tubular frame compatible with virtually any seating




Locking flip back arms with easy to release mechanism



ROHO® and Comfort available on-chair



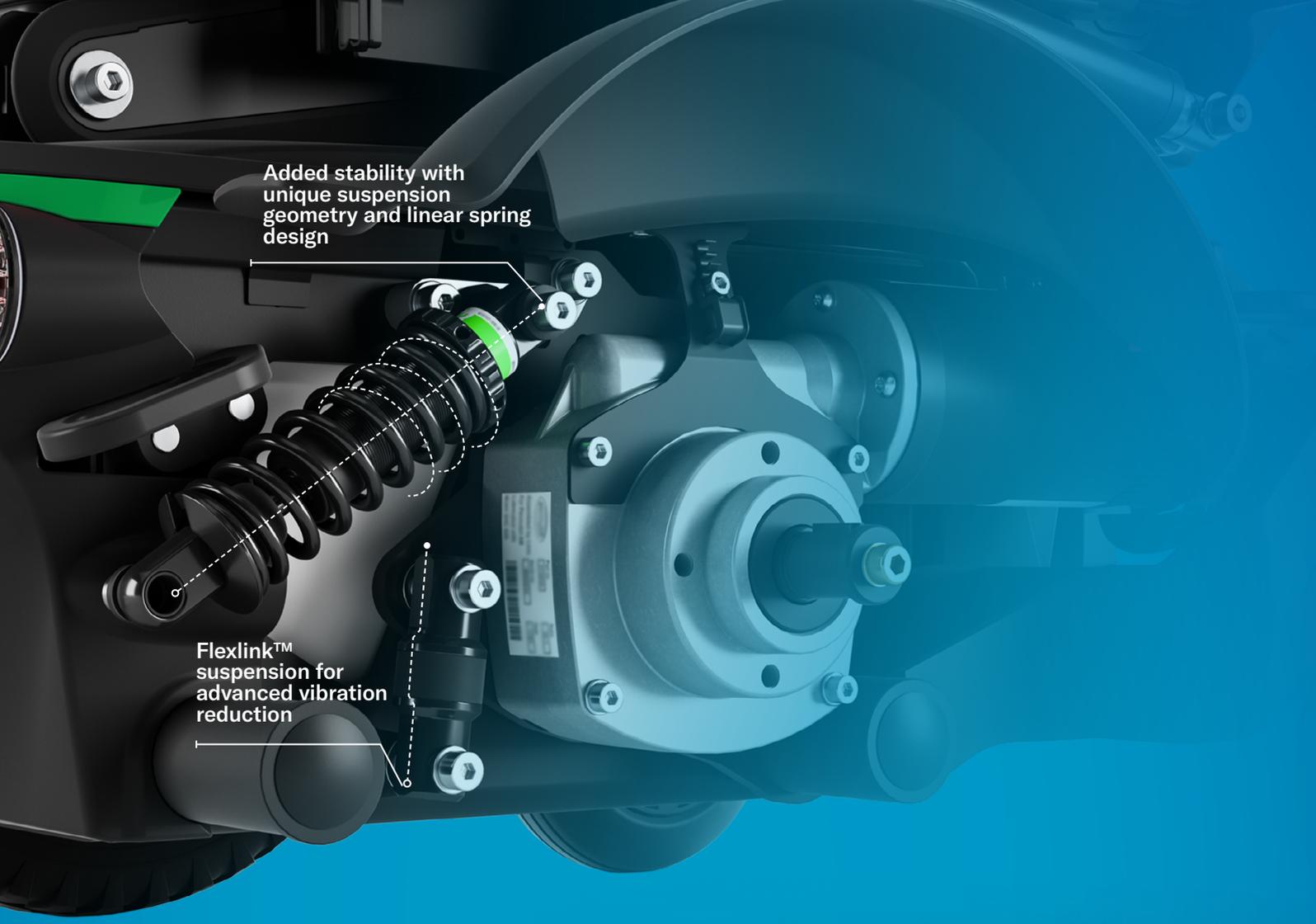

Added stability with
unique suspension
geometry and linear spring
design

Flexlink™
suspension for
advanced vibration
reduction

Proven technology

Built on the modular chassis of the Permobil M3 Corpus®, the M1 features the same reliable drive units, patented FlexLink™ suspension and precision engineered CasterGlide™ technology to help ensure smooth, confident operation in a variety of driving scenarios.



7" performance castor
design and adjustment
-free flutter dampening



Ready to go

Transport friendly features including standard crash tested transport brackets*, optional EZ fold backrest and a lower overall chair weight of 145 kg with batteries, expands available travel options.

*Crash Tested to ISO 7176-19:2008(E)

Specifications

Maximum user weight	136 kg ^{††}	Backrest height	18" – 29" (by 1" incr.)
Maximum speed	8 km/h	Seat depth	14" – 20" (by 1" incr.)
Range	25 km*	Seat width	14" – 20" (by 2" incr.)
Base width	24"	Armpad options	4"x13"
Base length	35.5"		
Base turning radius	20.5"		
Ground clearance	3"		
Recommended battery type	Group 34		
Weight including batteries	145 kg**		
Drive electronics	VR2 70A (Advanced)		
Seat-to-floor height	16"		
Tilt options	0° – 50° posterior		
Recline options	85° – 110° (manual)		
Legrest angle	80° – 135° (manual)		

†Specifications may vary based upon configuration. Please refer to your owners manual for warnings and instructions for the safe operation of your wheelchair in a variety of driving conditions.

††The maximum user weight for a seat depth of 14" cannot exceed 75 kg.

** Actual driving range will vary based on driving and battery conditions.*

***Includes weight of base, curved backrest and batteries.*

Colour Options



Sunburst Orange

Radiant Red

Ocean Blue

Power Pink

Galactic Green

Midnight Black