




The New Standard in Homecare



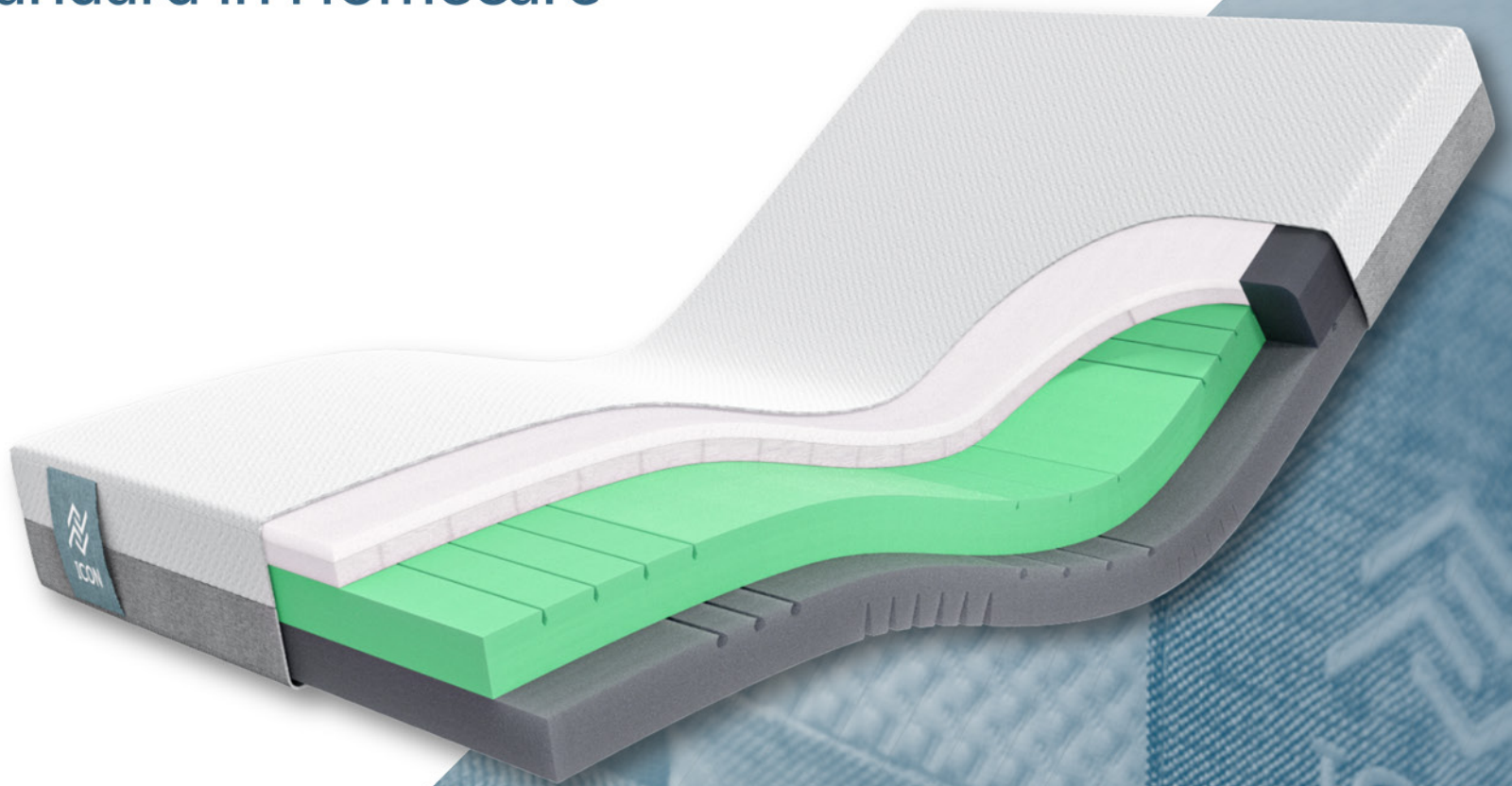

Designed for Safety, Engineered for Comfort

Universal Fit for Adjustable Beds

- Universal fitment – designed to fit and contour with adjustable bed frame
- 15 and 20cm standard thickness options with consistent comfort

Breathable & Comfortable for Homecare

- Integrated airflow channels for temperature and moisture control
- Medical-grade support with a home-like feel

Easier Repositioning & Increased Independence

- High support design to assist repositioning
- Soft, Medium or Firm options to support different comfort requirements

Safe & Supportive Transfers

- Firm, stabilising sides for safer transfers
- Radius edge for increased comfort and stability



Pressure Care Risk Rating

Shear & Friction	12345
Repositioning	12345
Infection Control	12345
Firmness	12345
Moisture & Microclimate	12345
Performance & Longevity	12345
Mobility & Function	12345
Pressure Care Rating	Low - Medium Risk
Warranty	5 Year Foam Core 2 Year Luxury Homecare Cover
Therapeutic Weight Loading	Soft (150kg) Medium (180kg) Firm (300kg)



Icon Mattress Features

Airflow & Temperature Regulation

- Open-cell Airoform layer for superior breathability.
- Side-to-side profiling promotes consistent ventilation.

Safety & Falls Prevention

- High-density strengthened sides to assist transfers.
- Reduces risk during transfers and repositioning.

Superior Repositioning Support

- Designed for users with reduced mobility.
- Triple Zone Pressure Care distributes pressure and stability where it's most needed.

Outstanding Comfort & Pressure Care

- Dual-layer Airoform foam for superior pressure care.
- Available in multiple comfort options for individual needs.



Foundational Support Foam

- Premium-grade firm foam enhances safety during transfers.
- Maintains support and comfort in varying user positions.

Advanced Medical Foam Technology

- 4-layer progressive support structure.
- Engineered for durability and long-term performance.

Adaptive Contouring for Adjustable Beds

- Unique profile hinge system for bed articulation.
- Channel profiling ensures a responsive and supportive feel.

Hygiene & Longevity

- Antimicrobial materials for enhanced infection control.
- High performance hospital-grade foams designed for long-term performance in a high-care environment.

Cover Options

The Icon range is available with a selection of high-performance covers designed for different care environments. Covers can also be purchased separately, allowing any mattress to be updated later to meet changing needs.

Luxury Quilted Homecare

- 4-way stretch, anti-shear fabric
- High breathability, assists cooling
- Homely quilted finish with medical-grade performance
- Supports effective pressure care
- Removable and replaceable
- Hard-wearing base built for adjustable bed movement



Advantiflex Classic

- 4-way stretch, anti-shear surface
- Waterproof, wipe-clean finish
- High-frequency sealed and welded seams
- Durable PVC base for added protection
- Compatible with disinfectants up to 1,000ppm




Premiflex Ultra

- 4-way stretch, anti-shear fabric for skin protection
- Waterproof, breathable surface, easy to wipe clean
- Engineered for demanding clinical environments
- Supports foam performance and hygiene management
- Compatible with disinfectants up to 10,000ppm




Safety, Function & Pressure Care

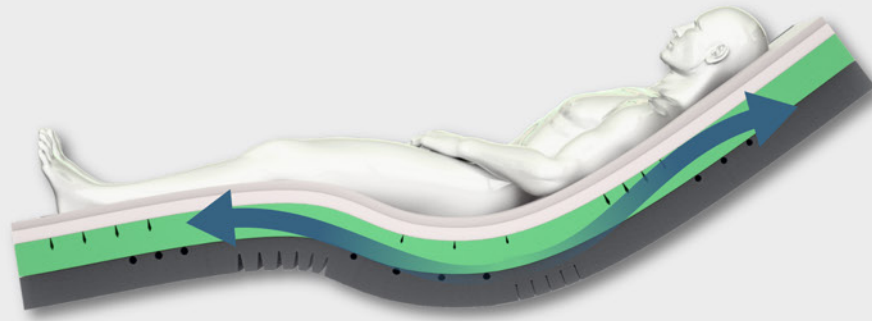

Breathability & Microclimate Control

- Side-to-side channels improve airflow and temperature regulation
- Open-cell foam layer increases ventilation
- Cover design supports microclimate management
- Responsive top balances comfort and pressure care

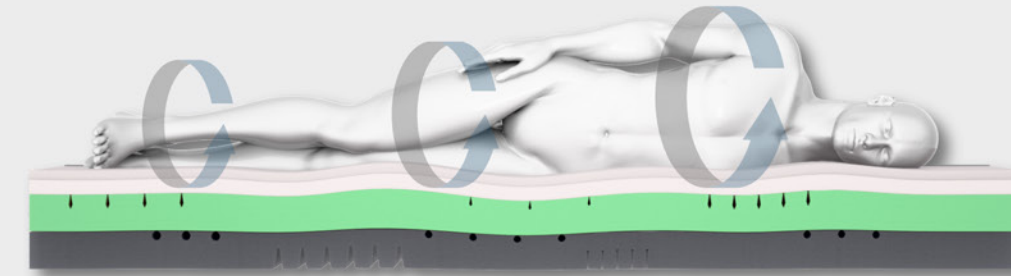

Ease of Repositioning

- High-grade comfort layers promote easier movement
- Airoform upper balances spring and immersion
- Reduced memory foam prevents sinking
- Triple Zone design supports natural pressure points


Safer Transfers & Stability

- Radius edge design promotes stability and confidence
- High-density sidewalls improve edge support
- Firm perimeter reduces risk of falls
- Enhanced stability supports safer transfers



Universal & Adaptive Fit

- Full hinge points maximise adjustable bed movement
- Knee and sacral profiling improve flexibility
- Compatible with all comfort options and thicknesses
- Comfort without compromising safety or durability

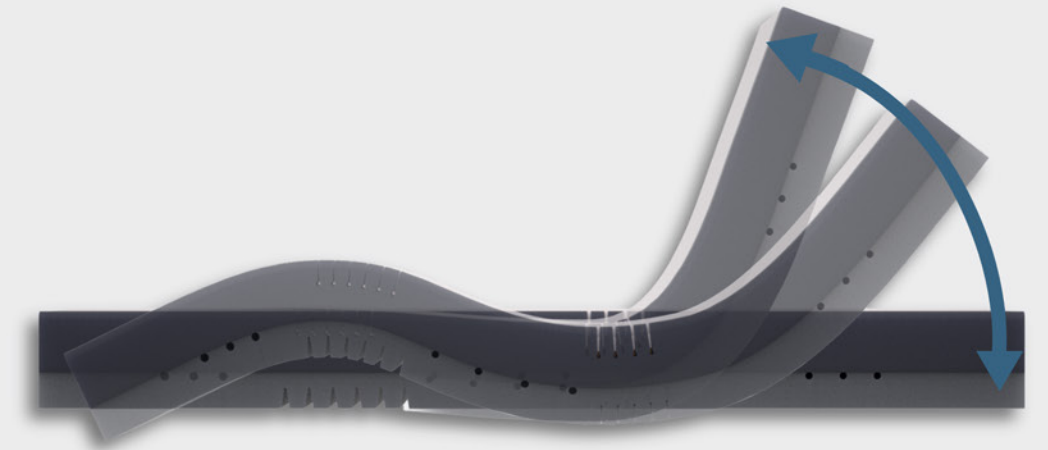

Comfort Options

The Icon mattress is available in Soft, Medium, and Firm comfort options, recognising that comfort is deeply personal and essential to effective clinical care. These options allow clinicians and caregivers to select a configuration that aligns with each client's unique needs, without compromising on the mattress's proven pressure care performance or core design principles.

ICON Soft

Ideal for clients who require gentle immersion and greater envelopment, helping reduce pressure peaks and enhance comfort for those with sensitivity or frailty.




ICON Medium

A balanced option that suits a broad range of clients, offering both responsive support and comfortable contouring.



ICON Firm


Designed for clients who need added stability or prefer a more supportive surface, promoting easier repositioning and transfer support.

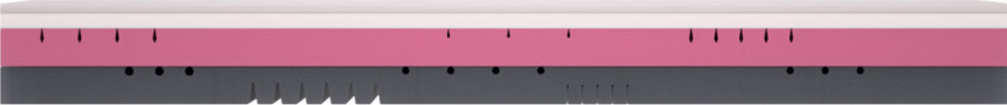

Height Options


The Icon mattress is offered in two height profiles—15cm and 20cm—to provide flexibility in both clinical functionality and user preference. Each option has been designed to support specific care needs while maintaining the Icon's high standards in pressure care and durability.

15cm Height

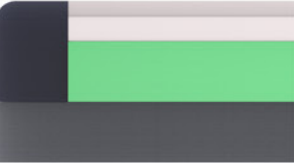

- Low-profile design helps reduce the risk of injury from bed-related falls
- Ideal for high-acuity settings and fall prevention protocols
- Seamless integration with profiling and adjustable beds
- Ensures correct alignment at key bed break points, such as the knee and hip
- Enables easier transfers for both clients and carers

20cm Height

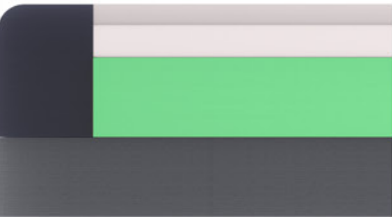

- Provides a more traditional mattress appearance and user experience
- Enhanced immersion and envelopment for improved pressure redistribution
- Supports higher weight capacities, making it suitable for bariatric care
- Increases long-term durability, ideal for high-use environments
- A versatile option for residential aged care and long-stay clinical settings



Stocked Icon Range

Other sizes and configurations available Made-to-Order

Comfort	Depth	Luxury Quilted	Advantiflex Cover
Soft	15 cm	Long Single, Long King Single	–
Soft	20 cm	Long Single, Long King Single, Long Double	Long King Single
Medium	15 cm	Long Single, Long King Single, Long Double	Long King Single, Long Double
Medium	20 cm	Single, Long Single, Long King Single, Long Split Queen	Long Single, Long King Single
Firm	15 cm	Long Single, Long King Single	–
Firm	20 cm	Long Single, Long King Single, Long Double, Long Queen	Long Single, Long King Single

Manufactured in Australia, Fully Customisable

Proudly crafted in Armidale, Australia, the Icon is available in all standard sizes and can also be tailored to your exact specifications, with adjustable length, width, and thickness, using Australian and European materials.

Scan the QR code to view the complete set of product codes.





The Mattress Makes The Bed™

Address

10 Ribbon Gum Ave, Armidale NSW 2350

Phone

1300 731 163

Email

sales@fortehealthcare.com.au

Web

fortehealthcare.com.au

Socials



© Copyright 2025 Forté Healthcare ICONV2 - 4/25

